



## PHARMACEUTICAL



# SP Technology


## Automation & Robotics



Multiple vials stored in preloaded trays are automatically loaded into our assembly system. Vials are weighed offline and then returned to the machine to be accurately filled using a servo driven pump. Devices are then transferred to a second weigh station to ensure the correct level of fluid has been dispensed. Auto capping is then followed by labelling on an SP Technology automatic labelling system for 360° label wrap, with full serialisation label inspection.

### Key features and equipment:

- Manual loading of prefilled trays (whilst machine running)
- Auto loading of devices from trays to machine
- Pre and post fill weighing
- Automated dispensing
- 360° Label wrapping & full serialisation
- Automated capping systems
- Output - 40 - 100 Parts per/minute

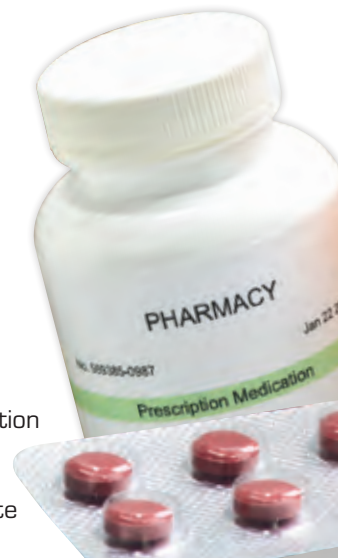

### SP Technology Limited

Camperdown Industrial Park, Dundee DD2 3SP  
Telephone: 01382 880088 | Email: [info@sptechology.co.uk](mailto:info@sptechology.co.uk)  
[www.sptechology.co.uk](http://www.sptechology.co.uk)

